



Specialty Devices

- 160 Disposal Basins
- 162 Drainage Bags
- 163 Fixation Device
- 164 Pen & Labels
- 165 Pressure Infusor Bag
- 166 Sharps Holders
- 167 Snares
- 169 Suture Retention
- 171 Transducers

Our specialty devices cover a wide range of quality products that provide unique solutions for your medical needs.

Disposal Basins



BackStop® Disposal Basin, 500ml fluid capacity and absorbent material. Includes particulate-free die cut sponge and seal-tight lid.

MBS100P



BackStop+™ Disposal Basin with temporary sharps holder, 500ml fluid capacity and absorbent material. Includes particulate-free die cut sponge and seal-tight lid.

MBSP100P



DugOut® Disposal Basin 1000ml fluid capacity and absorbent material. Includes particulate-free die cut sponge, secondary compartment for storage of ancillary products or saline and seal-tight lid.

DUG100P



DugOut® Disposal Basin with 1000ml fluid capacity and absorbent material. Includes particulate-free die cut sponge, secondary compartment for storage of ancillary products or saline and seal-tight lid with ShortStop® sharps holder.

DUG1SSP

Disposal Basins



MiniStop+™ Disposal Basin with temporary sharps holder, 100ml fluid capacity and absorbent material. Includes particulate-free die cut sponge and seal-tight lid.

MMSP100P



MiniStop® Disposal Basin, 100ml fluid capacity and absorbent material. Includes particulate-free die cut sponge and seal-tight lid.

MMS100P

Drainage Bags



Merit Drainage Depot™ 600ml Nephrostomy Bag, 24" (61cm). Kink-resistant adjustable tubing with fixed male luer, twist drain valve, VELCRO® brand fastener and additional male luer.

MDD600P




600ml bag, soft cloth backing, 24" (61cm) kink-resistant adjustable tubing with fixed male luer, twist drain valve and VELCRO® brand fastener. Additional male luer.

FZ624P

VELCRO® is a registered trademark of Velcro Industries B.V.

Fixation Device

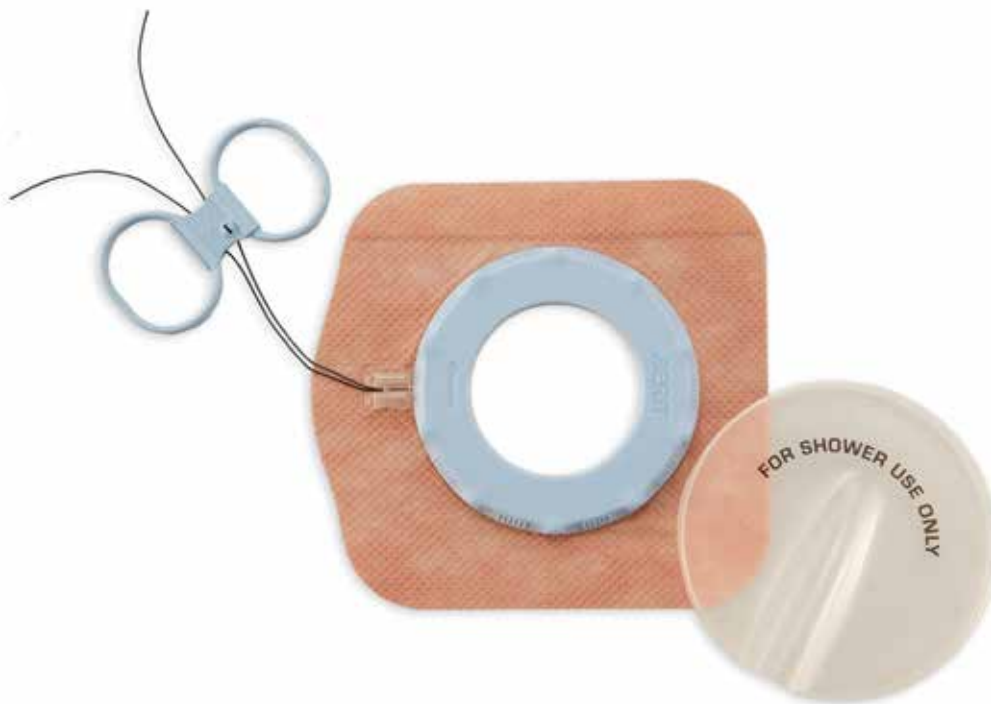

Revolution™ Needle-less Non-vascular
Catheter Securement Device (6F - 24F)

CFG624



Catalog Number	Description
PAL100	Contains pen and label set of 35
PAL200	Contains pen and label set of 24
PAL225	Contains pen, label set of 24 and guide wire bowl label set of 5
PAL250	Contains pen, label set of 24 and large bowl label set of 7
PAL300	Contains pen and label set of 12
PAL325	Contains pen, label set of 12 and guide wire bowl label set of 5
PAL350	Contains pen, label set of 12 and large bowl label set of 7
PAL400	Contains pen and blank label set of 12
PAL425	Contains pen, blank label set of 12 and guide wire bowl label set of 5
PAL450	Contains pen, blank label set of 12 and large bowl label set of 7
PAL500	Contains pen and guide wire bowl label set of 5
PAL600	Contains pen and large bowl label set of 7
PAL700	Contains pen, label set of 35 and large bowl label set of 7
K13-00476	Pen only

Pressure Infusor Bags

Many tubing sets available with non-DEHP material.

Color: White Material: Nylon Bag



500ml Pressure Infusor Bag

PIB500



1000ml Pressure Infusor Bag

PIB1000

SPECIALTY
DEVICES

Sharps Holders



ShortStop® Temporary
Sharps Holder

SS10P




ShortStop Advantage™ Temporary
Sharps Holder with Side Access Point

SS20P



ShortStop Grandstand® Temporary
Sharps Holder with Upright Access

SS30P

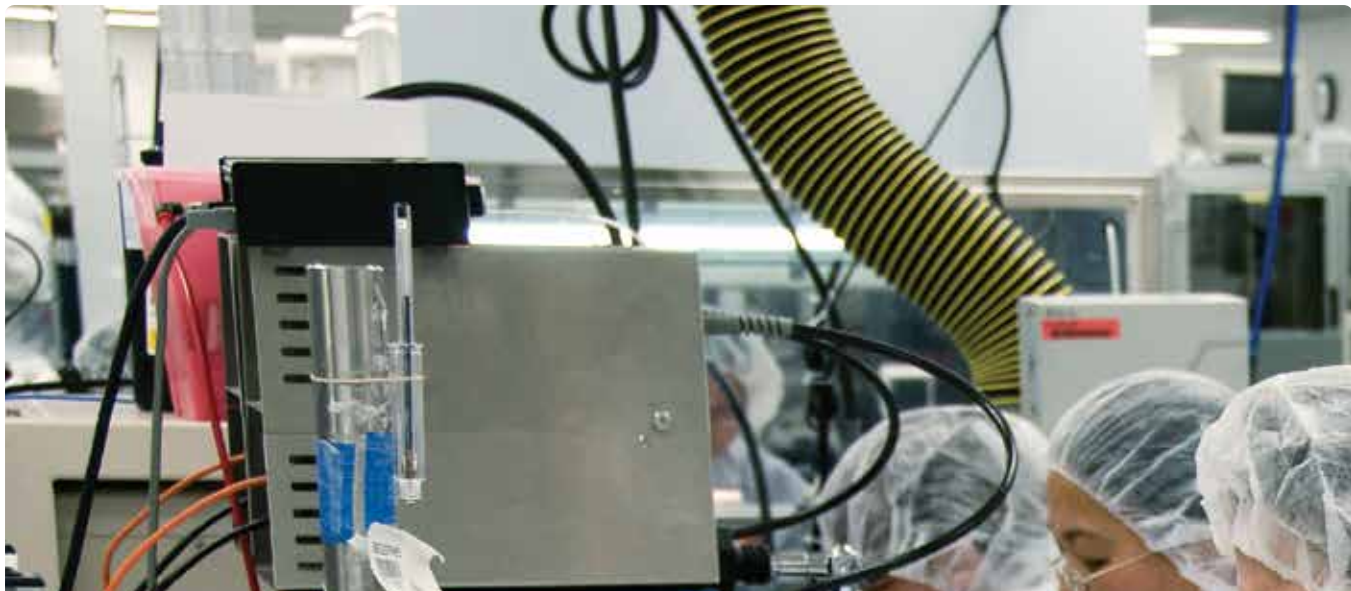

Snare

The EN Snare® Endovascular System is a multi-loop snare intended for use in the cardiovascular system or hollow viscus to retrieve and manipulate foreign objects.

The ONE Snare® Endovascular System is a single-loop snare designed to provide both accuracy and reliability for effective foreign body manipulation and retrieval.

SEE RELATED
PRODUCTS



Needles
145



Torque Devices
97



Hemostasis
Valves
243

May not be available for all OEM applications.
Contact a Merit OEM representative to discuss availability.

EN Snare® Endovascular System

Sold individually and includes a snare, torque device, introducer and snare catheter. Snare and catheter packaged separately.



The EN Snare® is designed to expand to vessel wall size for retrieval of hard to get side-walled objects.




The EN Snare® is designed with three loops to increase the probability of foreign body capture.



Interlaced Nitinol loops provide structural integrity when capturing objects.

Catalog Number	Working Diameter	Snare Length	French Size	Snare Catheter Length
EN1003004 Mini	2 - 4mm	175cm	3.2F	150cm
EN1003008 Mini	4 - 8mm	175cm	3.2F	150cm
EN2006010 Standard	6 - 10mm	120cm	6F	100cm
EN2006015 Standard	9 - 15mm	120cm	6F	100cm
EN2006020 Standard	12 - 20mm	120cm	6F	100cm
EN2007030 Standard	18 - 30mm	120cm	7F	100cm
EN2007045 Standard	27 - 45mm	120cm	7F	100cm

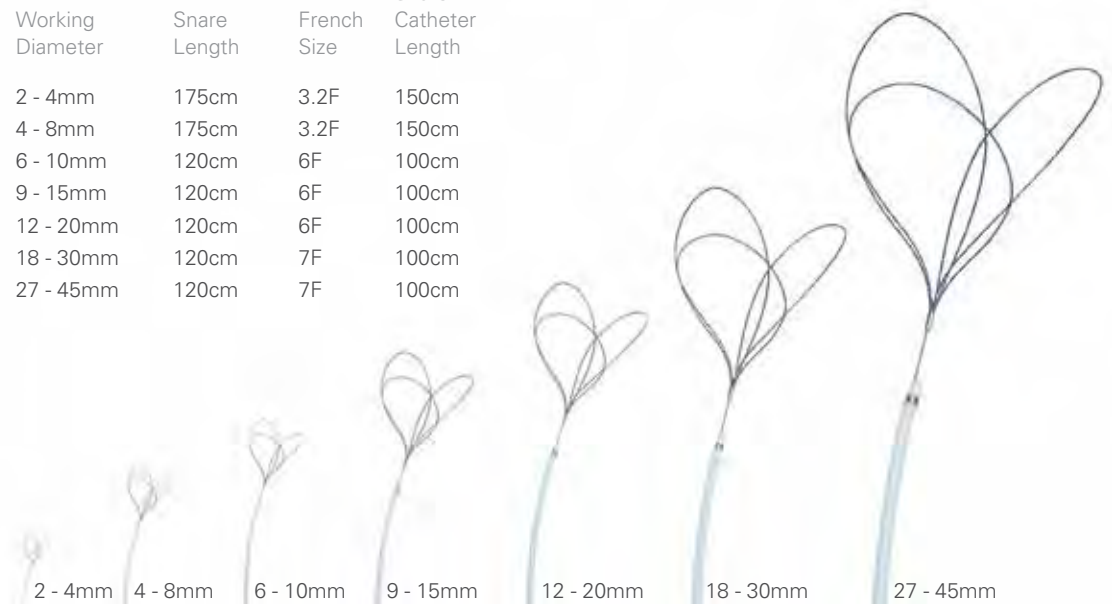

ONE Snare® Endovascular System

Catalog Number	Working Diameter	Snare Length	French Size	Snare Catheter Length
ONE500	5mm	120cm	4F	100cm
ONE1000	10mm	120cm	4F	100cm
ONE1500	15mm	120cm	6F	100cm
ONE2000	20mm	120cm	6F	100cm
ONE2500	25mm	120cm	6F	100cm
ONE3000	30mm	120cm	6F	100cm
ONE3500	35mm	120cm	6F	100cm
Petite Snare Kits				
ONE1001	10mm	65cm	4F	48cm
ONE2501	25mm	65cm	6F	48cm



SPECIALTY DEVICES



Suture Retention Device

A unique and effective method to secure a purse-string suture. Suture tension adjustments at the push of a button. Clear visibility of exit site. Enhanced patient comfort. Easy suture removal.

May not be available for all OEM applications.

Single use. Sold packaged sterile. For non-sterile and other packaging options, please contact a Merit OEM representative.

SEE RELATED PRODUCTS




Needles
145

Suture Retention Device


Easy, Step-by-Step Instructions



Place purse-string suture around the sheath; feed free ends of the suture through the snare loop of the device.

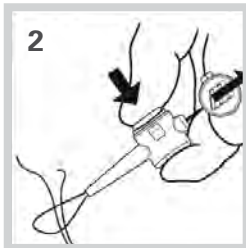

Depress the white button; pull the green tab to thread the free ends of the suture through the device. Discard green pull tab.

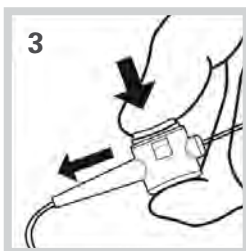

Depress the white button and slide the SLIP-NOT® forward to produce tension. Release the white button to secure the suture within the device.



SLIP-NOT®
SRD100

Transducers

- Reliable Precision
- Clear Fluid Path
- Customized Configurations also Available

SPECIFICATIONS

Operating Pressure Range	-50 to +300 mmHg
Sensitivity	5 μ V/V/mmHg, \pm 1%
Zero Offset	0 \pm 25 mmHg

Single use. Sold packaged non-sterile. For sterile and other packaging options, please contact a Merit OEM representative.

SEE RELATED
PRODUCTS



Syringes
195



Stopcocks
253



Manifolds
249



MERIT SENSOR

Merit Medical Systems' wholly owned subsidiary, Merit Sensor, designs and manufactures pressure sensors included in the Merit Transducer as well as Merit's digital inflation devices. See page 157 or go to www.MeritSensor.com

Meritrans® Blood Pressure Transducers



With Non-Winged Housing and 3-Way and 1-Way Stopcocks with 48" (122cm) cable

K11-00202P



With Non-Winged Housing and 3-Way Stopcock (Med Pressure, Rotating Male Adapter) and 1-Way Stopcock with 48" (122cm) cable

K11-00392AP



With Winged Housing and 2 Bonded Standard Pressure 1-Way Stopcocks with 12" (30cm) cable

K11-00098AP



With Non-Winged Housing and 3-Way Stopcock bonded to both ends with 48" (122cm) cable

K11-00144AP



With Non-Winged Housing and Bonded Standard Pressure 1-Way Stopcock with 48" (122cm) cable

K11-00137AP



With Winged Housing and Bonded Standard Pressure 3-Way and 1-Way Stopcocks with 12" (30cm) cable

K11-00100AP

SPECIALTY
DEVICES

Meritrans® Blood Pressure Transducers



With Non-Winged Housing with
12" (30cm) cable

MER212P



With Non-Winged Housing with
48" (122cm) cable

MER200P



With Winged Housing with
12" (30cm) cable

MER100P



With Winged Housing and 2 Bonded
Standard Pressure 3-Way Stopcocks with
12" (30cm) cable

K11-00099AP



With Winged Housing and Bonded
Standard Pressure 1-Way and 3-Way
Stopcocks with 12" (30cm) cable

K11-00096AP



With Winged Housing with
48" (122cm) cable

MER148P

SPECIALTY
DEVICES

Meritrans® Blood Pressure Transducers



Medium Pressure Transducer
0 - 100 PSI
501669001



High Pressure Transducer
0 - 440 PSI
501750001

Transducers Flush Device



Adult 3cc/hr flow rate.
White
270428001



Pediatric 30cc/hr flow rate.
Yellow
270428002