

Cultura™


TRANSPORT VIAL

Optimal Collection. Convenient Transport.

Introducing the new Cultura™ Transport Vial from Merit Medical. Developed as a critical part of streamlining the process of collecting and transporting specimen samples for a variety of clinical and laboratory applications. The features included in the Cultura™ Transport Vial help reduce the risk of specimen contamination and biohazard spills. Manufactured at Merit's Global Headquarters and Injection Molding Facility of Excellence in South Jordan, Utah—you can be confident in the stable supply of the highest quality product.

**MOLDED
IN USA**

Manufactured at Merit's corporate headquarters in South Jordan, Utah

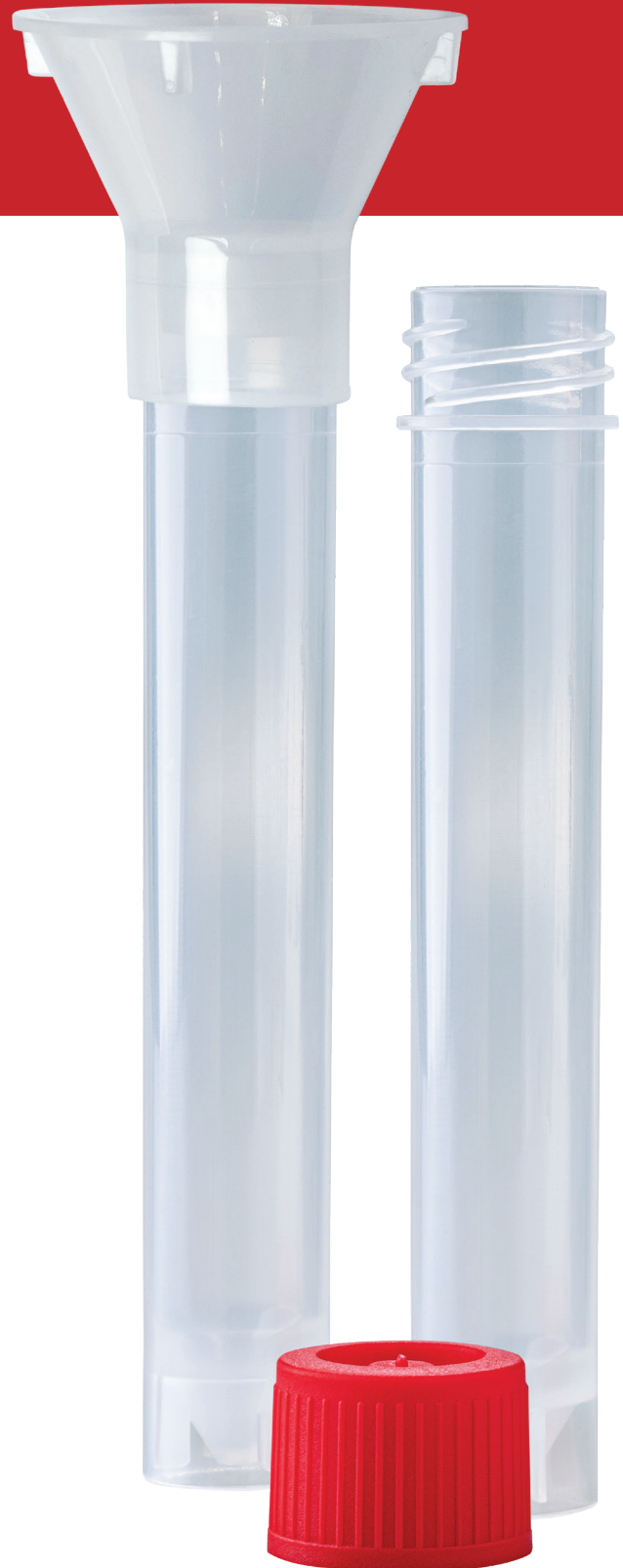

Cultura™

TRANSPORT VIAL

FEATURES & BENEFITS

- Self-standing flat bottom conical tube reduces the risk of biohazard spills
- Capture-cap design automatically grips swab as the vial is being sealed, reducing the risk of contamination when removing the swab
- Twist on cap with dual sealing surface design prevents leakage
- Saliva collection funnel attaches to vial with ease
- Available as a bulk non-sterile component or private labeled sterile kit
- Free of DNase, RNase, Human DNA, Pyrogens, and Endotoxins

SPECIFICATIONS

- Tube: Polypropylene, Cap: High Density Polyethylene, Funnel: Polypropylene
- Volume: 12mL
- Capture Cap: Holds 100 mm Swab
- Tube Dimensions: 16 mm x 101 mm

